


Screw Removing Instrument Kit (SRT10) Information And Instructions

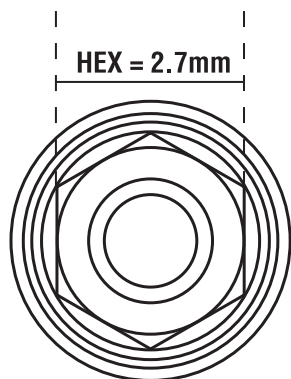

Hex –
Internal screw size 1.2mm

Product Description

The **COMET 3i** Screw Removing Instrument is designed to facilitate the removal of broken screws from dental implants with the additional benefit of not damaging the internal threads of the implant. The Screw Removing Instrument is packaged in a kit (SRT10).

Contents

The Screw Removing Instrument Kit (SRT10) Contains:

- (1) Drill Guide and Handle
- (1) 1.00mm diameter Reverse Latch Drill to fit a handpiece
- (1) 1.04mm diameter Left Hand Drill – Silver Handle
- (1) 1.18mm diameter Left Hand Drill – Bronze Handle
- (1) Hand Tap – Gold Handle

Indications

The Screw Removing Instrument Kit is used to retrieve screws that are broken below the occlusal surface of the implant and cannot be removed using other techniques. The instruments would be used when a burr, metal shaving or other defect on the broken surface of the screw prevents it from being easily removed.

Contraindications

There are no absolute contraindications however, broken screws, which protrude above the occlusal surface of the implant that can be gripped with a hemostat or slotted instrument, should be removed with such instruments. A broken screw is a salvage situation. The SRT10 or any other instrument may not be 100% effective in retrieving.

Caution: The drills and taps used in this kit are made of carbon steel. If a drill or tap is broken off inside the implant without being retrieved, removal of the implant may be necessary. Cold sterilization and careful drying will maintain instruments.

Directions For Use

1. If not already removed, the abutment or any other restorative components must be removed from the implant containing the broken screw. There must be an opening in the tissue of adequate size and shape to allow seating of a 4.5mm diameter abutment on the implant. The Drill Guide used is also 4.5mm in diameter.
2. Place the reverse drill with a latch type shank into a low speed handpiece capable of being regulated to approximately 200-400rpm. Make sure the handpiece motor is in reverse so that the drill rotates in a counter clockwise direction.
3. The handle end Drill Guide has an internal hex, which will fit over the male hex on the implant. Holding the handle in one hand, place the hexed end of the Drill Guide over the implant. With a light pressure on the handle, rotate the Drill Guide until it fits and is seated over the hex on the implant. Make sure the Drill Guide is held so that its longitudinal axis is parallel to the implant.



Screw Removing Instrument Kit
SRT10


Illustrations are representative and not to scale.

Screw Removing Instrument Kit (SRT10) Information And Instructions

(continued)

4. Holding the Drill Guide firmly over the implant and using copious external irrigation, insert the reverse drill through the hole on the superior end of the Drill Guide. With the handpiece rotating in reverse at 200-400rpm, continue to insert the drill until it engages the broken screw. Apply light to medium pressure on the handpiece so that the drill will begin to cut into the superior surface of the broken screw. Continue to apply irrigation during this step.
5. To avoid heat build-up in the implant; do not drill more than 8 to 10 seconds at a time, applying irrigation between drilling intervals to keep the implant cool. This process may have to be repeated several times to prepare a hole in the screw.
6. In many cases, the cutting action of the drill will cause the screw to move and begin to thread itself out of the implant. If this happens, reduce pressure on the handle and Drill Guide and apply only light pressure on the handpiece, allowing the screw to thread itself completely out of the implant.
7. Should the screw remain lodged in the implant, repeat the drilling described in the above paragraphs. The object is to drill a hole in the center of the broken screw approximately 1.0mm deep. The drill guide has an internal stop, which prevents over penetration of the drill.
8. If the screw remains in place after a hole has been drilled through the longitudinal axis of the broken screw, the Hand Drills are designed to provide additional torque for final removal. Use the 1.04mm Hand Drill (silver colored handle) first. If it will not grip the screw due to the internal hole being too large, substitute the 1.18mm Hand Drill (bronze colored handle).
9. Insert the drill end into the internal implant thread, being careful to engage the hole previously drilled in the broken screw. Apply light pressure in the apical direction and slowly rotate the Hand Drill counter clockwise. The drill should wedge itself into the previously generated hole and apply additional removal torque to the screw, removing the screw from the implant. Please note that the success of this step depends upon the drill entering the hole in the screw. If the 1.04mm Drill fails to remove the screw and it appears to be too small to wedge into the drilled hole, the 1.18mm Drill may be used.
10. Should both of the hand drills fail to remove the screw, other instruments may be tried. For example, a small round burr slightly larger than the drilled hole may be used successfully.
11. Once the screw has been removed, it is important to examine the internal threads of the implant to determine if these have been damaged. The best method of performing this inspection is to thread an Implant Guide Pin (WSK30) into the implant. The Guide Pin should thread into the implant with little or no resistance realizing that some resistance may occur because of minor damage to the threads.
12. If significant thread damage has occurred, a tap is included in the SRT10 Kit for retapping the implant. The tap must be used very carefully since cross-threading it into the implant would permanently damage the threads, which could result in the implant having to be removed.




This instrument is not guaranteed to be capable of removing all broken abutment screws and should be used with extreme care. The decision to use this instrument is entirely at the discretion of the practitioner. For technical support or more information please contact  at 800-342-5454.

BIOMET 3i LLC. makes no other warranty, express or implied, except that its products shall be free from defects in material and/or workmanship and shall be of merchantable quality. This warranty applies only to the original purchaser. In the event of a defect, please notify BIOMET 3i of the defect prior to returning this instrument, BIOMET 3i will, at its option, either repair, replace or issue credit for such defective product. Purchaser assumes all risks and liability resulting from the use of these products, whether used separately or in combination with other products not manufactured by BIOMET 3i.

Providing Solutions – One Patient At A Time is a registered trademark of BIOMET 3i LLC. BIOMET 3i and design are trademarks and BIOMET is a registered trademark of BIOMET, Inc. ©2016 BIOMET 3i LLC. All rights reserved.

CE
0086

P-IIS032
REV B 05/16

 BIOMET 3i
4555 Riverside Drive
Palm Beach Gardens, FL 33410
1.800.342.5454
Outside U.S. + 561.776.6700
Fax: + 561.776.1272
www.biomet3i.com

 BIOMET 3i Dental Iberica S.L.
WTC Almeda Park, Ed. 4, Planta 2
C/Tirso de Molina, 40
08940 - Cornellà de Llobregat
(Barcelona) Spain
Phone: + 34 934 705 500
Fax: + 34 933 717 849


PROVIDING SOLUTIONS – ONE PATIENT AT A TIME™